

Ligand Exchange Chromatography Columns

Lists summarizing elution volumes of various saccharides using Shodex columns is available. Please refer to our website (<http://www.shodex.com/en/>) or technical notebook (No.2 and 3).

Features

- SC1011**
 - Separates saccharides by combination of ligand exchange and size exclusion modes
 - SC1821**
 - Three types of counter ions are available: Ca²⁺, Pb²⁺, and Na⁺
 - SP0810**
 - Only water is required for the analysis of neutral sugars
 - KS-801 to 802**
 - SC1011 and SC1821 fulfill USP L19 and L22 requirements
 - SP0810 fulfills USP L22 and L34 requirements
 - KS-801 and KS-802 fulfill USP L22 and L58 requirements
-
- KS-803 to 806**
 - Suitable for separation of polysaccharides by size exclusion mode
 - Can be used in combination with other columns e.g., KS-802 and KS-801
 - Only water is required for the analysis of neutral sugars
 - Fulfills USP L22 and L58 requirements
-
- DC-613**
 - Separates elements by combination of ligand exchange and HILIC modes
 - SZ5532**
 - DC-613 can analyze sugars without removing sodium salts in the sample
 - SC1211**
 - SZ5532 is recommended for the separation of disaccharides or trisaccharides
 - SC1211 is suitable for separating sugar alcohols
 - DC-613 fulfills USP L22 and L58 requirements
 - SZ5532 fulfills USP L22 requirements
 - SC1211 fulfills USP L19 and L22 requirements
-
- SC1011-7F**
 - Fulfills mannitol analysis requirements of JP, USP, and EP methods
 - Ca²⁺ modified ligand exchange chromatography column
 - Only water is required for the analysis of neutral sugars
 - Fulfills USP L19 and L22 requirements

Standard columns

[Ligand exchange and size exclusion]

Product Code	Product Name	Plate Number (TP/column)	Functional Group (Counter Ion)	Exclusion Limit (Pullulan)	Particle Size (μm)	Column Size (mm) I.D. x Length	Shipping Solvent
F6378102	SUGAR SC1011	≥ 13,000	Sulfo (Ca ²⁺)	1,000	6	8.0 x 300	H ₂ O
F6378103	SUGAR SC1821	≥ 13,000	Sulfo (Ca ²⁺)	10,000	6	8.0 x 300	H ₂ O
F6700090	SUGAR SC-G 6B (SUGAR SC-LG)	(guard column)	Sulfo (Ca ²⁺)	–	10	6.0 x 50	H ₂ O
F6378105	SUGAR SP0810	≥ 11,000	Sulfo (Pb ²⁺)	1,000	7	8.0 x 300	H ₂ O
F6700081	SUGAR SP-G 6B (SUGAR SP-G)	(guard column)	Sulfo (Pb ²⁺)	–	10	6.0 x 50	H ₂ O
F6378010	SUGAR KS-801	≥ 17,000	Sulfo (Na ⁺)	1,000	6	8.0 x 300	H ₂ O
F6378020	SUGAR KS-802	≥ 17,000	Sulfo (Na ⁺)	10,000	6	8.0 x 300	H ₂ O
F6378025	SUGAR KS-803	≥ 17,000	Sulfo (Na ⁺)	50,000	6	8.0 x 300	H ₂ O
F6378035	SUGAR KS-804	≥ 17,000	Sulfo (Na ⁺)	400,000	7	8.0 x 300	H ₂ O
F6378050	SUGAR KS-805	≥ 9,000	Sulfo (Na ⁺)	5,000,000	17	8.0 x 300	H ₂ O
F6378060	SUGAR KS-806	≥ 9,000	Sulfo (Na ⁺)	*(50,000,000)	17	8.0 x 300	H ₂ O
F6700020	SUGAR KS-G 6B (SUGAR KS-G)	(guard column)	Sulfo (Na ⁺)	–	10	6.0 x 50	H ₂ O

*() Estimated value Base Material: Styrene divinylbenzene copolymer

[Ligand exchange and HILIC]

Product Code	Product Name	Plate Number (TP/column)	Functional Group (Counter Ion)	Particle Size (μm)	Pore Size (Å)	Column Size (mm) I.D. x Length	Shipping Solvent
F7001003	RSpak DC-613	≥ 5,500	Sulfo (Na ⁺)	6	100	6.0 x 150	H ₂ O/CH ₃ CN=30/70
F6700170	RSpak DC-G 4A (RSpak DC-G)	(guard column)	Sulfo (Na ⁺)	10	–	4.6 x 10	H ₂ O/CH ₃ CN=30/70
F7001300	SUGAR SZ5532	≥ 5,500	Sulfo (Zn ²⁺)	6	–	6.0 x 150	H ₂ O/CH ₃ CN=30/70
F6700110	SUGAR SZ-G	(guard column)	Sulfo (Zn ²⁺)	6	–	4.6 x 10	H ₂ O/CH ₃ CN=30/70
F7001400	SUGAR SC1211	≥ 5,500	Sulfo (Ca ²⁺)	6	50	6.0 x 250	H ₂ O/CH ₃ CN=75/25
F6700120	SUGAR SC1211G 4A (SUGAR SC-G)	(guard column)	Sulfo (Ca ²⁺)	10	–	4.6 x 10	H ₂ O/CH ₃ CN=75/25

Base Material: Styrene divinylbenzene copolymer

■ For mannitol analysis following JP, USP, and EP methods

Product Code	Product Name	Functional Group (Counter Ion)	Particle Size (μm)	Column Size (mm) I.D. x Length	Shipping Solvent
F6379300	EP SC1011-7F	Sulfo (Ca ²⁺)	8	7.8 × 300	H ₂ O
F6700090	SUGAR SC-G 6B (SUGAR SC-LG) (guard column)	Sulfo (Ca ²⁺)	10	6.0 × 50	H ₂ O
F6379230	USPpak MN-431	Sulfo (Ca ²⁺)	8	4.0 × 250	H ₂ O

*See page 79 for USP40-NF35 Column List.

Base Material: Styrene divinylbenzene copolymer

■ Preparative columns *Preparative columns are made to order.

[Ligand exchange and size exclusion]

Product Code	Product Name	Plate Number (TP/column)	Particle Size (μm)	Column Size (mm) I.D. x Length	Standard Column
F6502007	SUGAR KS-2001	≥ 7,000	13	20.0 × 300	KS-801
F6502008	SUGAR KS-2002	≥ 7,000	13	20.0 × 300	KS-802
F6502009	SUGAR KS-2003	≥ 8,000	13	20.0 × 300	KS-803
F6502010	SUGAR KS-2004	≥ 6,000	18	20.0 × 300	KS-804
F6502011	SUGAR KS-2005	≥ 6,000	18	20.0 × 300	KS-805
F6502012	SUGAR KS-2006	≥ 6,000	18	20.0 × 300	KS-806
F6700002	SUGAR KS-G 8B (SUGAR KS-LG)	(guard column)	13	8.0 × 50	(guard column)

[Ligand exchange and HILIC]

Product Code	Product Name	Plate Number (TP/column)	Particle Size (μm)	Column Size (mm) I.D. x Length	Standard Column
F6514013	RSPak DC-2013	≥ 6,000	10	20.0 × 300	DC-613
F6700402	RSPak DC-G 8B (RSPak DC-LG)	(guard column)	10	8.0 × 50	(guard column)

Elution volumes of saccharides analyzed by Shodex columns

[Partial list only; refer to our website for complete list]

Substances	Elution Volume (mL)					
	SP0810	SC1011	KS-801	SZ5532	NH2P-50 4E	SC1211
Arabinose	10.42	8.91	8.21	5.11	6.18	5.56
D-Arabitol	15.86	11.33	7.63	7.27	6.29	8.16
Dulcitol	20.18	12.76	7.40	9.46	7.45	11.28
meso-Erythritol	12.70	10.09	7.86	5.73	5.43	6.27
D(-)-Fructose	11.05	8.85	7.71	5.37	6.75	5.90
D(+)-Fucose	10.48	8.84	8.09	4.50	5.43	4.96
D(+)-Galactose	9.74	7.98	7.58	6.46	8.10	4.98
Gentiobiose	7.22	6.08	5.75	10.50	16.36	*
Glucose	8.63	7.30	7.17	5.87	8.61	4.76
myo-Inositol	12.77	8.86	7.99	12.63	9.96	7.87
Isomaltose	7.68	6.26	5.95	10.57	15.18	*
Isomaltotriose	7.09	5.75	5.34	21.17	27.55	*
1-Kestose	6.79	5.75	5.26	13.09	20.11	*
Kojibiose	7.56	6.21	5.88	9.65	14.82	*
Lactitol	13.27	8.09	6.13	16.35	11.82	6.67
Lactose	8.05	6.51	5.99	10.12	13.27	4.07
Lactulose	9.13	6.99	6.19	9.16	10.72	4.65
Maltitol	12.23	8.26	6.03	13.04	11.82	6.77
Maltose	7.85	6.34	5.94	8.67	14.24	*
Maltotriose	7.48	5.89	5.38	13.79	24.96	*
Mannitol	15.80	11.10	7.23	8.75	7.39	9.03

Substances	Elution Volume (mL)					
	SP0810	SC1011	KS-801	SZ5532	NH2P-50 4E	SC1211
D-Mannose	10.72	8.17	7.64	5.83	7.84	5.01
Melibiose	8.16	6.45	5.98	11.69	14.70	4.23
Nystose	6.38	5.45	4.93	20.05	31.90	*
Palatinin	2peaks	2peaks	5.90	2peaks	12.73	2peaks
Palatinose	7.84	6.45	5.89	8.08	12.12	3.99
Panose	7.14	5.78	5.32	16.87	25.60	*
D(+)-Raffinose	7.14	5.78	5.29	16.36	20.25	*
Rhamnose	9.77	8.23	7.37	3.93	5.52	4.43
D(-)-Ribose	19.35	13.66	9.04	4.82	5.45	8.64
D(-)-Sorbitol	21.61	13.31	7.42	9.79	7.09	11.88
Sorbose	9.67	8.03	7.38	5.12	7.35	4.92
Stachyose	6.82	5.57	4.97	—	36.22	*
Sucrose	7.54	6.29	5.87	7.91	11.87	*
α-D-Talose	21.33	12.59	8.76	5.69	6.47	8.51
Trehalose	7.62	6.27	5.78	10.85	13.25	*
Trehalulose	8.92	6.95	6.10	9.54	11.68	4.78
Xylitol	19.87	13.14	7.94	7.77	6.10	10.16
Xylobiose	8.16	6.68	6.40	5.65	9.05	*
D(+)-Xylose	9.21	7.90	7.71	4.55	6.58	4.48
D-Xylose	10.64	9.02	8.04	4.06	5.41	5.07

(-) → Not detected (-) → Overlap with solvent peak

(-) → Not detected (-) → Overlap with solvent peak

Column : SUGAR SP0810, SC1011, KS-801
 Eluent : H₂O
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 80°C


Column : SUGAR SC1211
 Eluent : H₂O/CH₃CN=65/35
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 70°C

Column : SUGAR SZ5532
 Eluent : H₂O/CH₃CN=25/75
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 60°C

Column : Asahipak NH2P-50 4E
 Eluent : H₂O/CH₃CN=25/75
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 30°C


Saccharides anomer separation

Saccharides may present their anomers at lower temperatures. By setting the SUGAR series columns at higher temperatures will prevent the anomer separation and this results in providing better chromatograms of each saccharide. Sample : 0.5% each, 10µL

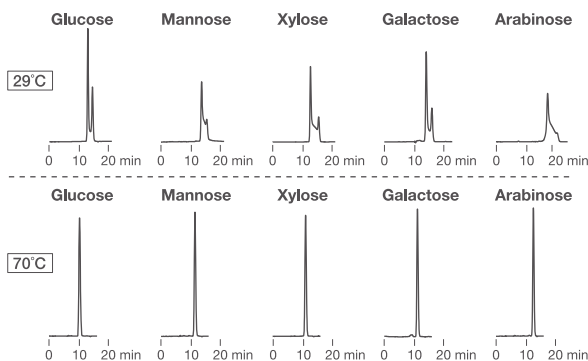


Column : Shodex SUGAR SC1011
 Eluent : H₂O
 Flow rate : 0.7mL/min
 Detector : RI
 Column temp.: 29°C, 70°C

Calibration curves for KS-800 series using pullulan

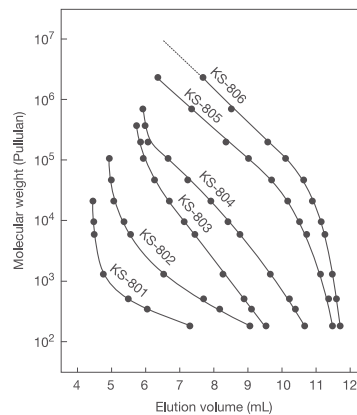


Column : Shodex SUGAR KS-800 series
 Eluent : H₂O
 Detector : RI
 Column temp. : 80°C

Hydrolyzed starch

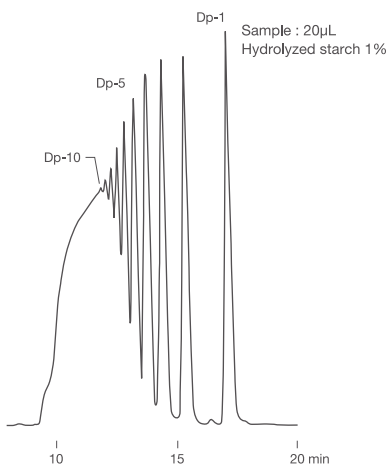

Sample : 20µL
 Hydrolyzed starch 1%
 Column : Shodex SUGAR KS-802 x 2
 Eluent : H₂O
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 80°C


Biomass sugars

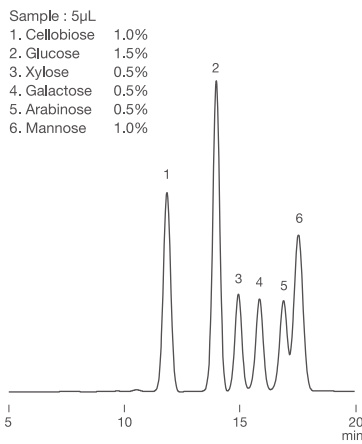

Sample : 5µL
 1. Cellobiose 1.0%
 2. Glucose 1.5%
 3. Xylose 0.5%
 4. Galactose 0.5%
 5. Arabinose 0.5%
 6. Mannose 1.0%

Column : Shodex SUGAR SP0810
 Eluent : H₂O
 Flow rate : 0.6mL/min
 Detector : RI
 Column temp. : 85°C


Ketohexoses

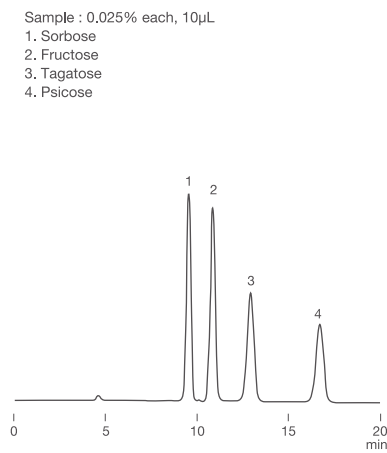

Sample : 0.025% each, 10µL
 1. Sorbose
 2. Fructose
 3. Tagatose
 4. Psicose

Column : Shodex SUGAR SP0810
 Eluent : H₂O
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 80°C


Pinitol

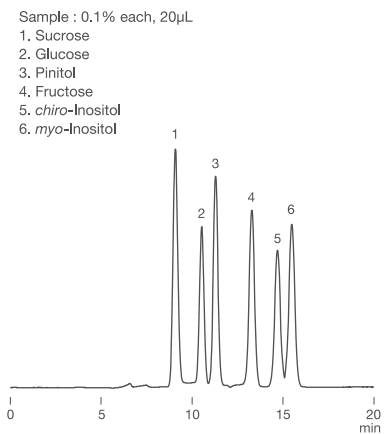

Sample : 0.1% each, 20µL
 1. Sucrose
 2. Glucose
 3. Pinitol
 4. Fructose
 5. chiro-Inositol
 6. myo-Inositol

Column : Shodex SUGAR SP0810
 Eluent : H₂O
 Flow rate : 0.8mL/min
 Detector : RI
 Column temp. : 85°C


Oligosaccharides in soybean

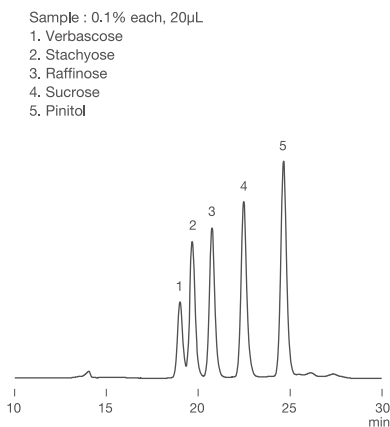

Sample : 0.1% each, 20µL
 1. Verbascose
 2. Stachyose
 3. Raffinose
 4. Sucrose
 5. Pinitol

Column : Shodex SUGAR KS-802 + KS-801
 Eluent : H₂O
 Flow rate : 0.6mL/min
 Detector : RI
 Column temp. : 85°C

Saccharides related to raffinose biosynthesis

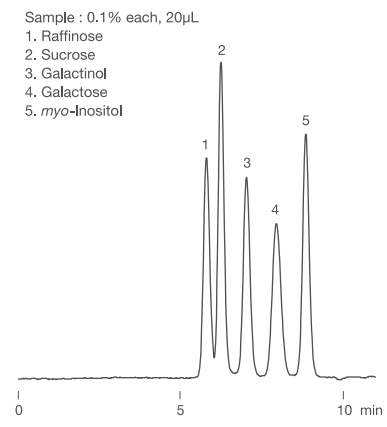

Sample : 0.1% each, 20µL
 1. Verbascose
 2. Sucrose
 3. Galactinol
 4. Galactose
 5. myo-Inositol

Column : Shodex SUGAR SC1011
 Eluent : H₂O
 Flow rate : 1.0mL/min
 Detector : RI
 Column temp. : 80°C

